

**Instrument Model:** 

Ion Source:

## **Instrument Tuning Report**

**Instrument Name:** ZenoTOF™ 7600 System LC200

ZenoTOF 7600 System

OptiFlow 1-50µL Micro/MicroCal

Serial Number: FB24962403

Manufacturer: AB Sciex

 Target Instrument:
 ZenoTOF 7600 System

Optimization Set: Positive Quick Status Check

## 1. Positive MS Check

**Channel Alignment** 


Step Result Details

Achieve Stable Spray / Modify Pass

Slope: 0.000490178399948099, previous was

0.000490178399948099

Delay: -13.4570174949306, previous was -13.4570174949306

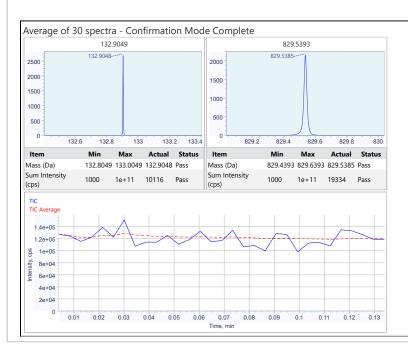

**Pass** 

## Sum of 10 spectra 132.9049 829.5393 132.9044-829.5381-2 1.5e+05 1e+05 1e+05 5e+04 132.8 133 133.2 133.4 829 2 829.4 829.6 Item Max Mass (Da) 132.6549 133.1549 132.9044 Mass (Da) 829.2893 829.7893 829.5381 Resolution 0 100000 38972 Pass Resolution 0 100000 43442 Pass Sum Intensity Sum Intensity 1e+11 6.4518e+05 Pass 1e+11 2.0225e+06 Pass (cps) Align to channel: 2 Result: Pass Min Max Initial Channel 1 vs. 2 Slope -50 50 -0.127427973 -0.118501576 Pass Channel 1 vs. 2 Intercept -50 50 -10.78771791 -11.453747983 Pass Channel 3 vs. 2 Slope 0.072660388 0.075212626 Pass -3.996284445 -4.025329752 Pass Channel 4 vs. 2 Slope 0.111825314 0.115082975 Channel 4 vs. 2 Intercept -50 50 9.059756629 8.616777451

www.SCIEX.com Page 1 of 2



**Instrument Model:** 

Ion Source:

Serial Number:

## **Instrument Tuning Report**

ZenoTOF™ 7600 System LC200

ZenoTOF 7600 System

OptiFlow 1-50µL Micro/MicroCal

FB24962403


Manufacturer: AB Sciex

ZenoTOF 7600 System **Target Instrument:** 

**Optimization Set:** Positive Quick Status Check

TOF MS Mass Check **Pass** Slope: 0.000490178385123102, previous was 0.000490178399948099

Delay: -13.4577071296833, previous was -13.4570174949306

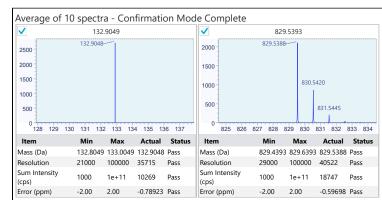

TOF MS/MS Mass Check **Pass** Slope: 0.000490177658387613, previous was 0.000490177859332247 Delay: -13.4789145092723, previous was -13.4646724934156 Average of 10 spectra - Confirmation Mode Complete 185.1285 298.2125 494.3337 607.4178 2000 1 5e+04 185.1285 494.3339 607.4178-298.2126 8000 1.5e+04 1500 6000 1e+04 1000 1e+04 4000 611.4127 500 5000 187,1440 498.3286 490 491 492 493 494 495 496 497 498 499 181 182 183 184 185 186 187 188 189 190 294 295 296 297 298 299 300 301 302 303 603 604 605 606 607 608 609 610 611 612 Item Min Max Actual Status Item Max Actual Status Item Min Max Actual Status Item Min Max Actual 185.0285 185.2285 185.1285 Pass 298.1125 298.3125 298.2126 Pass 607.3178 607.5178 607.4178 Mass (Da) Mass (Da) Mass (Da) 494.2337 494.4337 494.3339 Pass Mass (Da) Pass Resolution 21000 100000 27713 Pass Resolution 100000 36261 Pass 25000 100000 42230 Resolution 26000 100000 40070 22000 Resolution Pass Pass Sum Intensity Sum Intensity Sum Intensity Sum Intensity 100 1e+11 11543 Pass 100 1e+11 48413 100 1e+11 1.1755e+05 Pass 100 1e+11 1.0782e+05 Pass (cps) (cps) (cps) -0.01156 Pass -2.00 2.00 -2.00 2.00 0.49313 Pass -2.00 2.00 0.45199 Pass -2.00 2.00 -0.05063 Pass Error (ppm) Error (ppm) Error (ppm) Error (ppm)

TOF MS/MS Mass Check Zeno On Slope: 0.000490179362685442, previous was **Pass** 0.000490178955931249 Delay: -13.4430372035668, previous was -13.4350160450402 Average of 10 spectra - Confirmation Mode Complete 185.1285 298.2125 494.3337 607.4178 85 1285 4000 3000 800 2500 2500 3000 600 2000 2000 1500 611.4129 400 2000 1500 1000 1000 200 1000 500 187.1439 295.2378 181 182 183 184 185 186 187 188 189 190 294 295 296 297 298 299 300 301 302 303 490 491 492 493 494 495 496 497 498 499 603 604 605 606 607 608 609 610 611 612 Max Item Max Actual Status Item Max Actual Status Min Actual Status Min Max Actual Status Min Min Item Item Mass (Da) 185 0285 185 2285 185 1285 Pass Mass (Da) 298 1125 298 3125 298 2126 Pass Mass (Da) 494 2337 494 4337 494 3340 Pass Mass (Da) 607 3178 607 5178 607 4179 Pass Resolution 21000 100000 33481 Pass Resolution 22000 100000 41051 Pass Resolution 25000 100000 38447 Pass Resolution 26000 100000 40905 Pass Sum Intensity Sum Intensity Sum Intensity Sum Intensity 1e+11 100 4244.8 Pass 100 13952 Pass 100 26882 100 21766 Pass (cps) (cps) (cps) (cps) Error (ppm) Error (ppm) Error (ppm) Error (ppm

www.SCIEX.com Page 2 of 2