

## **Instrument Tuning Report**

**nstrument Name:** ZenoTOF™ 7600 System LC200

Instrument Model: ZenoTOF 7600 System

**Ion Source:** OptiFlow 1-50µL Micro/MicroCal

0.005 0.01 0.015 0.02 0.025 0.03 0.035 0.04 0.045 0.05 0.055 0.06 0.065 0.07 0.075 0.08 0.085 0.09 0.095 0.1


Serial Number: FB24962403


Manufacturer: AB Sciex

Target Instrument:ZenoTOF 7600 System

**Optimization Set:** Positive Quick Status Check

## 1. Positive MS Check

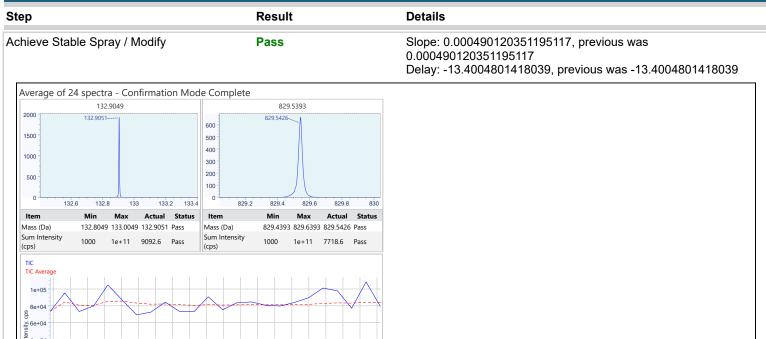

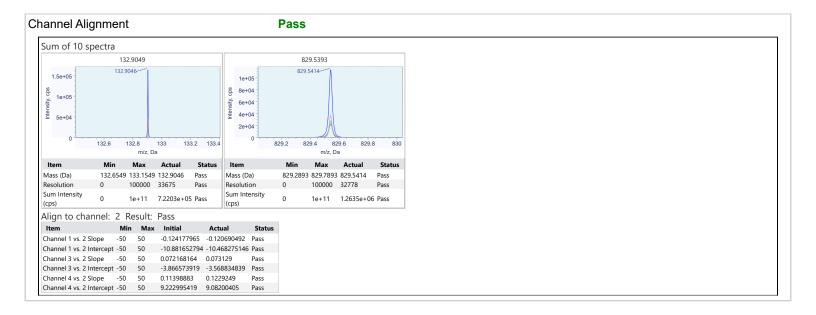

www.SCIEX.com Page 1 of 2



**Instrument Model:** 

Error (ppm)

-2.00

2.00

0.26278 Pass

## **Instrument Tuning Report**

-2.00

2.00

Error (ppm)

ZenoTOF™ 7600 System LC200

ZenoTOF 7600 System

ZenoTOF 7600 System **Target Instrument:** 

Manufacturer:

OptiFlow 1-50µL Micro/MicroCal Ion Source: **Optimization Set:** Positive Quick Status Check

Serial Number: FB24962403

TOF MS Mass Check **Pass** Slope: 0.000490119109922009, previous was 0.000490120351195117

Delay: -13.4435457615753, previous was -13.4004801418039

(cps)


Error (ppm)

-2.00

2.00

0.33499 Pass

AB Sciex

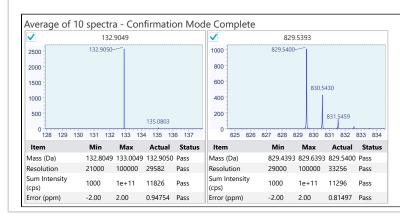

TOF MS/MS Mass Check **Pass** Slope: 0.000490117958121145, previous was 0.000490119979250965 Delay: -13.4769930435253, previous was -13.4471429657627 Average of 10 spectra - Confirmation Mode Complete 298.2125 494.3337 185 1285 298 2128 494 3341-607 4180 5000 800 3000 6000 2500 4000 600 2000 4000 3000 1500 400 2000 1000 200 498.3287 1000 500 295.2380 294 295 296 297 298 299 300 301 302 303 181 182 183 184 185 186 187 188 189 190 490 491 492 493 494 495 496 497 498 499 603 604 605 606 607 608 609 610 611 612 Max Max Actual Item Min Actual Item Min Max Item Item 185.0285 185.2285 185.1285 Pass 494.2337 494.4337 494.3341 Pass 607.3178 607.5178 607.4180 Pass Mass (Da) Mass (Da) 298.1125 298.3125 298.2128 Pass Mass (Da) Mass (Da) Resolution 100000 29646 Pass 21000 100000 23486 Pass Resolution 22000 100000 26932 Pass Resolution 25000 Resolution 26000 100000 28369 Pass Sum Intensity Sum Intensity Sum Intensity Sum Intensity 100 1e+11 6328.3 Pass 100 1e+11 26357 Pass 100 1e+11 64519 Pass 100 1e+11 58769 Pass

Error (ppm)

-2.00

2.00

0.87686 Pass

1.00067 Pass

TOF MS/MS Mass Check Zeno On **Pass** Slope: 0.000490119171212364, previous was 0.000490120837774385 Delay: -13.455220836718, previous was -13.4565806511591 Average of 10 spectra - Confirmation Mode Complete 185.1285 298.2125 494.3337 607.4178 298.2124~ 494.3338-400 2000 1000 300 1500 1000 200 611.4124 1000 500 100 500 498.3284 183 184 186 187 188 189 190 295 296 297 298 299 300 301 302 303 492 493 494 495 498 499 605 606 607 608 609 610 611 612 185 Item Min Max Actual Status Item Max Actual Status Item Min Max Min Max Actual Status Mass (Da) 185.0285 185.2285 185.1284 Pass Mass (Da) 298.1125 298.3125 298.2124 Pass Mass (Da) 494.2337 494.4337 494.3338 Pass Mass (Da) 607.3178 607.5178 607.4176 Pass Resolution 21000 100000 25752 Resolution 100000 30901 Pass Resolution 25000 100000 31900 Resolution 26000 100000 30997 Pass Pass Sum Intensity Sum Intensity Sum Intensity Sum Intensity 14864 (cps) (cps) (cps) -2.00 2.00 -0.4049 Pass -2.00 2.00 -0.48319 Pass -2.00 2.00 0.11076 Pass -2.00 2.00 -0.33594 Pass Error (ppm) Error (ppm) Error (ppm) Error (ppm)

www.SCIEX.com Page 2 of 2