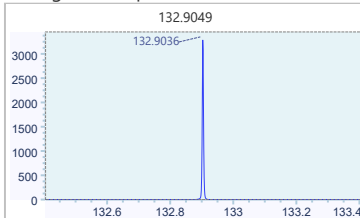
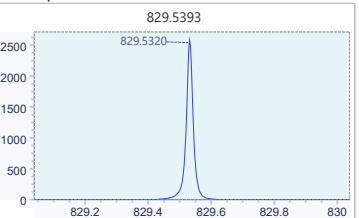
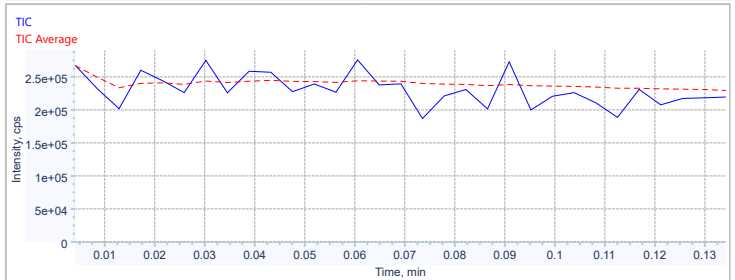


Instrument Name: ZenoTOF™ 7600 System
Instrument Model: ZenoTOF 7600 System
Ion Source: OptiFlow 50-200µL Micro/MicroCal
Serial Number: FB24962403

Manufacturer: AB Sciex
Target Instrument: ZenoTOF 7600 System
Optimization Set: Positive Quick Status Check

1. Positive MS Check

Step	Result	Details																														
Achieve Stable Spray / Modify	Pass	Slope: 0.00049010248950908, previous was 0.00049010248950908 Delay: -13.4517952035879, previous was -13.4517952035879																														
Average of 30 spectra - Confirmation Mode Complete																																
																																
<table><tr><th>Item</th><th>Min</th><th>Max</th><th>Actual</th><th>Status</th></tr><tr><td>Mass (Da)</td><td>132.8049</td><td>133.0049</td><td>132.9036</td><td>Pass</td></tr><tr><td>Sum Intensity (cps)</td><td>1000</td><td>1e+11</td><td>14560</td><td>Pass</td></tr></table>	Item	Min	Max	Actual	Status	Mass (Da)	132.8049	133.0049	132.9036	Pass	Sum Intensity (cps)	1000	1e+11	14560	Pass	<table><tr><th>Item</th><th>Min</th><th>Max</th><th>Actual</th><th>Status</th></tr><tr><td>Mass (Da)</td><td>829.4393</td><td>829.6393</td><td>829.5320</td><td>Pass</td></tr><tr><td>Sum Intensity (cps)</td><td>1000</td><td>1e+11</td><td>24638</td><td>Pass</td></tr></table>		Item	Min	Max	Actual	Status	Mass (Da)	829.4393	829.6393	829.5320	Pass	Sum Intensity (cps)	1000	1e+11	24638	Pass
Item	Min	Max	Actual	Status																												
Mass (Da)	132.8049	133.0049	132.9036	Pass																												
Sum Intensity (cps)	1000	1e+11	14560	Pass																												
Item	Min	Max	Actual	Status																												
Mass (Da)	829.4393	829.6393	829.5320	Pass																												
Sum Intensity (cps)	1000	1e+11	24638	Pass																												
																																

Channel Alignment

Pass

Sum of 10 spectra

132.9049

Intensity, cps

2e+05

1.5e+05

1e+05

5e+04

0

132.6

132.8

133

133.2

133.4

m/z, Da

132.9031

Item	Min	Max	Actual	Status
Mass (Da)	132.6549	133.1549	132.9031	Pass
Resolution	0	100000	36220	Pass
Sum Intensity (cps)	0	1e+11	7.8859e+05	Pass

829.5393

Intensity, cps

2e+05

1.5e+05

1e+05

5e+04

0

829.2

829.4

829.6

829.8

830

m/z, Da

829.5306

Item	Min	Max	Actual	Status
Mass (Da)	829.2893	829.7893	829.5306	Pass
Resolution	0	100000	41294	Pass
Sum Intensity (cps)	0	1e+11	2.0239e+06	Pass

Align to channel: 2 Result: Pass

Item	Min	Max	Initial	Actual	Status
Channel 1 vs. 2 Slope	-50	50	-0.125456582	-0.122518284	Pass
Channel 1 vs. 2 Intercept	-50	50	-10.753156757	-10.531321103	Pass
Channel 3 vs. 2 Slope	-50	50	0.073699958	0.076863231	Pass
Channel 3 vs. 2 Intercept	-50	50	-3.749719316	-4.02318819	Pass
Channel 4 vs. 2 Slope	-50	50	0.1195739	0.126990006	Pass
Channel 4 vs. 2 Intercept	-50	50	8.702500278	8.814334773	Pass

Instrument Name: ZenoTOF™ 7600 System
Instrument Model: ZenoTOF 7600 System
Ion Source: OptiFlow 50-200µL Micro/MicroCal
Serial Number: FB24962403


Manufacturer: AB Sciex
Target Instrument: ZenoTOF 7600 System
Optimization Set: Positive Quick Status Check

TOF MS Mass Check

Pass

Slope: 0.000490104374824947, previous was 0.00049010248950908
Delay: -13.462872515078, previous was -13.4517952035879

Average of 10 spectra - Confirmation Mode Complete

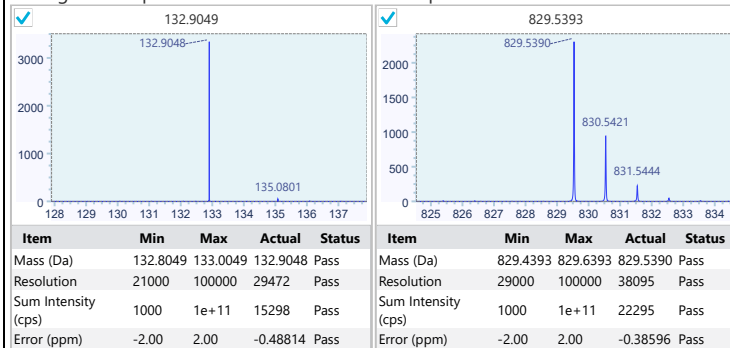


TOF MS/MS Mass Check

Pass

Slope: 0.000490103455115775, previous was 0.00049010551950563
Delay: -13.4920110822169, previous was -13.4808324288948

Average of 10 spectra - Confirmation Mode Complete

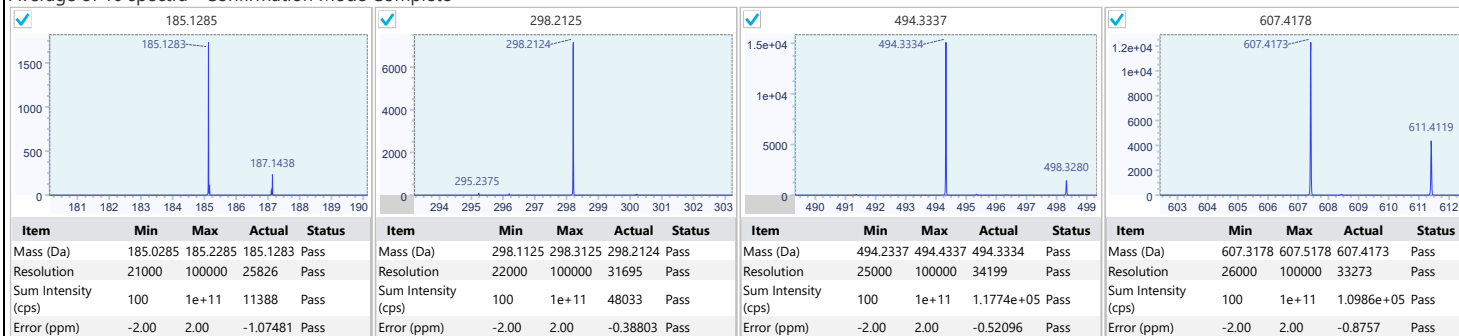

TOF MS/MS Mass Check Zeno On

Pass

Slope: 0.000490104958223297, previous was 0.000490102775985005
Delay: -13.4656955555172, previous was -13.464094087671

Average of 10 spectra - Confirmation Mode Complete

