

Instrument Tuning Report

nstrument Name: ZenoTOF™ 7600 System

Instrument Model: ZenoTOF 7600 System

Ion Source: OptiFlow 50-200µL Micro/MicroCal

Serial Number: FB24962403

Manufacturer: AB Sciex

 Target Instrument:
 ZenoTOF 7600 System

Optimization Set: Positive Quick Status Check

1. Positive MS Check


Step Result Details

Achieve Stable Spray / Modify

Pass Slope: 0.000490104715917106, previous was

0.000490104715917106

Delay: -13.4520306749997, previous was -13.4520306749997

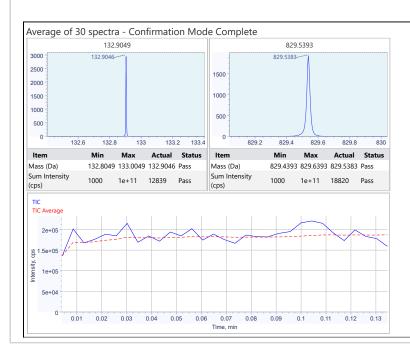

Pass Channel Alignment Sum of 10 spectra 132.9049 829.5393 132.9042-829.5372 නු 1.5e+05 1.5e+05 1e+05 1e+05 132.8 133 133.2 133.4 829 2 829.4 829.6 Item Max Mass (Da) 132.6549 133.1549 132.9042 Mass (Da) 829.2893 829.7893 829.5372 Resolution 0 100000 35606 Pass Resolution 0 100000 39278 Sum Intensity Sum Intensity 1e+11 7.6197e+05 Pass 0 1e+11 1.8144e+06 Pass (cps) Align to channel: 2 Result: Pass Min Max Initial Actual Channel 1 vs. 2 Slope -50 50 -0.123566735 -0.12521812 Pass Channel 1 vs. 2 Intercept -50 50 -10.680500117 -10.257904335 Pass Channel 3 vs. 2 Slope 0.073143452 0.075747727 Pass -3.555585587 -3.836817511 Pass Channel 4 vs. 2 Slope 0.125310911 0.128812663 Channel 4 vs. 2 Intercept -50 50 9 110417281 8 764401363

www.SCIEX.com Page 1 of 2



Instrument Model:

Ion Source:

Instrument Tuning Report

ZenoTOF™ 7600 System

ZenoTOF 7600 System

Target Instrument: OptiFlow 50-200µL Micro/MicroCal

Optimization Set: Positive Quick Status Check

AB Sciex

Serial Number: FB24962403

TOF MS Mass Check **Pass**


Slope: 0.000490104999979427, previous was

0.000490104715917106

Manufacturer:

Delay: -13.4656499865628, previous was -13.4520306749997


ZenoTOF 7600 System

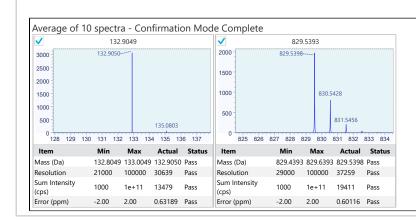

TOF MS/MS Mass Check **Pass** Slope: 0.000490103899074133, previous was

0.000490103485882046


Delay: -13.489058160333, previous was -13.4865775742939

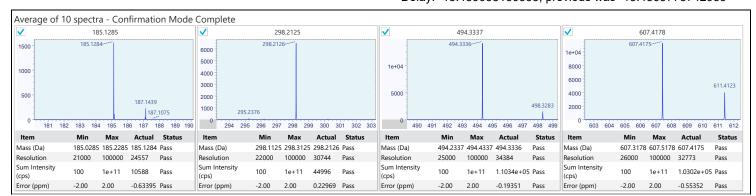

TOF MS/MS Mass Check Zeno On Slope: 0.000490105275456648, previous was **Pass**

0.000490104800406062

Delay: -13.4721739798397, previous was -13.4645892335978

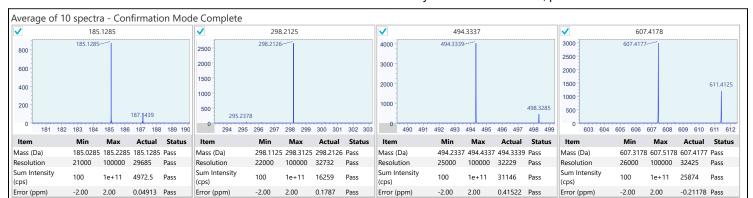

www.SCIEX.com Page 2 of 2