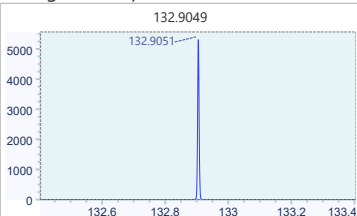
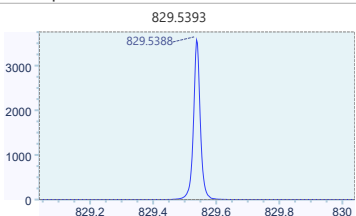
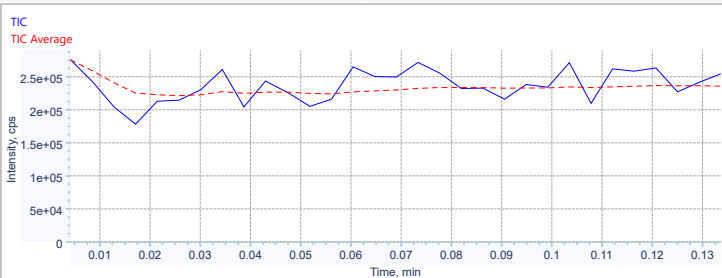


Instrument Name: ZenoTOF™ 7600 System
Instrument Model: ZenoTOF 7600 System
Ion Source: OptiFlow 50-200µL Micro/MicroCal
Serial Number: FB24962403

Manufacturer: AB Sciex
Target Instrument: ZenoTOF 7600 System
Optimization Set: Positive Quick Status Check

1. Positive MS Check

Step	Result	Details																														
Achieve Stable Spray / Modify	Pass	Slope: 0.000490174949761718, previous was 0.000490174949761718 Delay: -13.4887403415621, previous was -13.4887403415621																														
Average of 30 spectra - Confirmation Mode Complete																																
																																
<table><tr><th>Item</th><th>Min</th><th>Max</th><th>Actual</th><th>Status</th></tr><tr><td>Mass (Da)</td><td>132.8049</td><td>133.0049</td><td>132.9051</td><td>Pass</td></tr><tr><td>Sum Intensity (cps)</td><td>1000</td><td>1e+11</td><td>26138</td><td>Pass</td></tr></table>	Item	Min	Max	Actual	Status	Mass (Da)	132.8049	133.0049	132.9051	Pass	Sum Intensity (cps)	1000	1e+11	26138	Pass	<table><tr><th>Item</th><th>Min</th><th>Max</th><th>Actual</th><th>Status</th></tr><tr><td>Mass (Da)</td><td>829.4393</td><td>829.6393</td><td>829.5388</td><td>Pass</td></tr><tr><td>Sum Intensity (cps)</td><td>1000</td><td>1e+11</td><td>32016</td><td>Pass</td></tr></table>	Item	Min	Max	Actual	Status	Mass (Da)	829.4393	829.6393	829.5388	Pass	Sum Intensity (cps)	1000	1e+11	32016	Pass	
Item	Min	Max	Actual	Status																												
Mass (Da)	132.8049	133.0049	132.9051	Pass																												
Sum Intensity (cps)	1000	1e+11	26138	Pass																												
Item	Min	Max	Actual	Status																												
Mass (Da)	829.4393	829.6393	829.5388	Pass																												
Sum Intensity (cps)	1000	1e+11	32016	Pass																												
																																

Channel Alignment

Pass

Sum of 10 spectra

Item	Min	Max	Actual	Status	Item	Min	Max	Actual	Status
Mass (Da)	132.6549	133.1549	132.9044	Pass	Mass (Da)	829.2893	829.7893	829.5379	Pass
Resolution	0	100000	41099	Pass	Resolution	0	100000	43147	Pass
Sum Intensity (cps)	0	1e+11	7.1471e+05	Pass	Sum Intensity (cps)	0	1e+11	2.0584e+06	Pass
Align to channel: 2 Result: Pass									
Item	Min	Max	Initial	Actual	Status				
Channel 1 vs. 2 Slope	-50	50	-0.050293963	-0.039423508	Pass				
Channel 1 vs. 2 Intercept	-50	50	-10.206346413	-10.99007258	Pass				
Channel 3 vs. 2 Slope	-50	50	-0.005097428	0.002874032	Pass				
Channel 3 vs. 2 Intercept	-50	50	-3.637824259	-4.900827005	Pass				
Channel 4 vs. 2 Slope	-50	50	-0.046731276	-0.03631523	Pass				
Channel 4 vs. 2 Intercept	-50	50	8.197220641	7.59231876	Pass				

Instrument Name: ZenoTOF™ 7600 System
Instrument Model: ZenoTOF 7600 System
Ion Source: OptiFlow 50-200µL Micro/MicroCal
Serial Number: FB24962403


Manufacturer: AB Sciex
Target Instrument: ZenoTOF 7600 System
Optimization Set: Positive Quick Status Check

TOF MS Mass Check

Pass

Slope: 0.000490175279530462, previous was 0.000490174949761718
Delay: -13.4711490518821, previous was -13.4887403415621

Average of 10 spectra - Confirmation Mode Complete

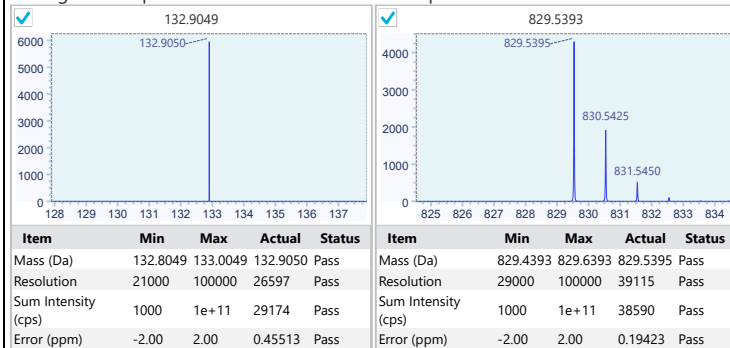


TOF MS/MS Mass Check

Pass

Slope: 0.000490174010308048, previous was 0.000490174386218911
Delay: -13.5046576161516, previous was -13.5002542809143

Average of 10 spectra - Confirmation Mode Complete

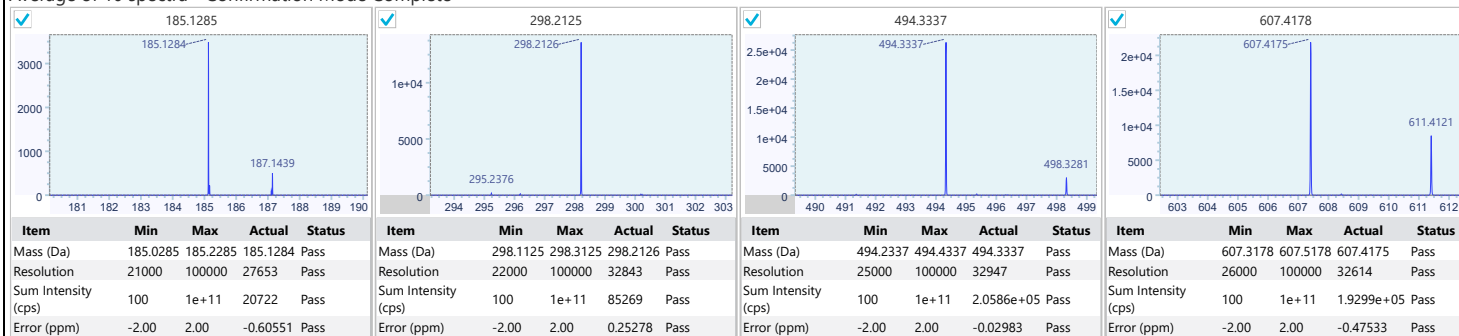

TOF MS/MS Mass Check Zeno On

Pass

Slope: 0.000490174489367797, previous was 0.000490174948048446
Delay: -13.4813644978154, previous was -13.4855328463341

Average of 10 spectra - Confirmation Mode Complete

