

Ion Source:

## **Instrument Tuning Report**

OptiFlow 1-50µL Micro/MicroCal

ZenoTOF™ 7600 System

ZenoTOF 7600 System **Instrument Model:** 

FB24962403 Serial Number:


Manufacturer: **AB Sciex** 

**Target Instrument:** ZenoTOF 7600 System

**Optimization Set:** Positive Quick Status Check

## 1. Positive MS Check

Step Result **Details** Achieve Stable Spray / Modify **Pass** Slope: 0.000490122440795605, previous was 0.000490122440795605 Delay: -13.4149201866171, previous was -13.4149201866171 Average of 30 spectra - Confirmation Mode Complete

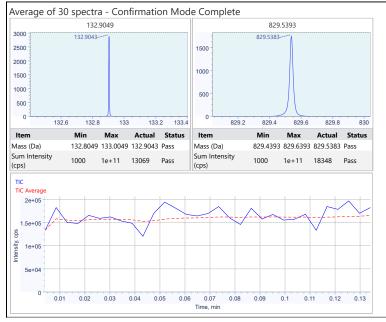

## **Pass Channel Alignment** Sum of 10 spectra 132.9049 829.5393 132.9036 829,5360 1.5e+05 1.5e+05 1e+05 5e+04 5e+04 132.8 133 133.2 133.4 829 2 829.4 829.6 Item Max Actual Mass (Da) 132.6549 133.1549 132.9036 Mass (Da) 829.2893 829.7893 829.536 Resolution 0 100000 34985 Pass Resolution 0 100000 37281 Pass Sum Intensity Sum Intensity 1e+11 7.4809e+05 Pass 0 1e+11 1.8587e+06 Pass (cps) Align to channel: 2 Result: Pass Min Max Initial Actual Channel 1 vs. 2 Slope -50 50 -0.120439599 -0.122429401 Pass Channel 1 vs. 2 Intercept -50 50 -11.657194339 -10.62482212 Pass Channel 3 vs. 2 Slope 0.073094164 0.075104962 -3.640378945 -3.878255432 Pass Channel 4 vs. 2 Slope 0.126218514 0.120371098 Channel 4 vs. 2 Intercept -50 50 8 212593487 8 803329125

www.SCIEX.com Page 1 of 2



**Instrument Model:** 

## **Instrument Tuning Report**

ZenoTOF™ 7600 System

ZenoTOF 7600 System

**Target Instrument:** 

ZenoTOF 7600 System

OptiFlow 1-50µL Micro/MicroCal Ion Source:

**Optimization Set:** Positive Quick Status Check

Serial Number: FB24962403


TOF MS Mass Check **Pass**  Slope: 0.000490122481511774, previous was

AB Sciex

0.000490122440795605

Manufacturer:


Delay: -13.4739355687954, previous was -13.4149201866171

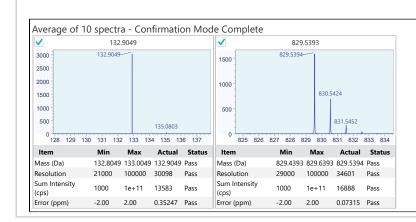

Slope: 0.000490121547319674, previous was TOF MS/MS Mass Check **Pass** 

0.00049012149313815


Delay: -13.4898539070873, previous was -13.4942891211385

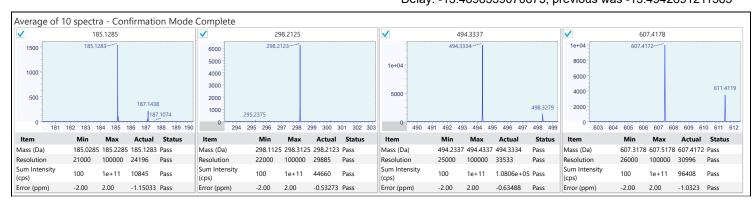

TOF MS/MS Mass Check Zeno On Slope: 0.000490122716558554, previous was **Pass** 

0.000490122965199601

Delay: -13.4967942816746, previous was -13.4769215801117

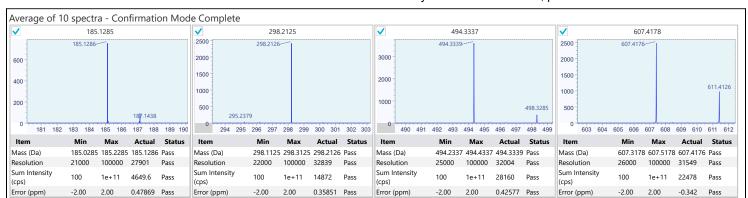

www.SCIEX.com Page 2 of 2