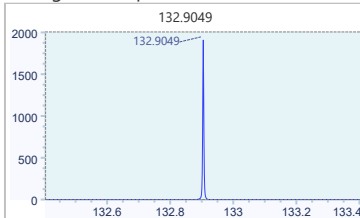
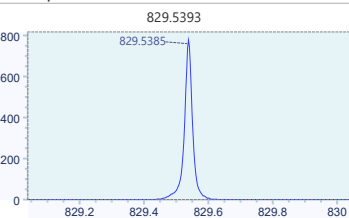
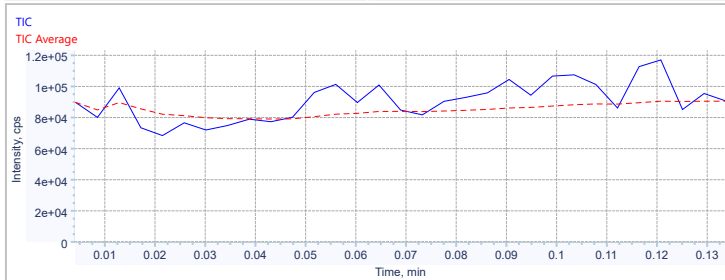


Instrument Name: ZenoTOF™ 7600 System
Instrument Model: ZenoTOF 7600 System
Ion Source: OptiFlow 50-200µL Micro/MicroCal
Serial Number: FB24962403

Manufacturer: AB Sciex
Target Instrument: ZenoTOF 7600 System
Optimization Set: Positive Quick Status Check



1. Positive MS Check

Step	Result	Details																														
Achieve Stable Spray / Modify	Pass	Slope: 0.000490188612503516, previous was 0.000490188612503516 Delay: -13.474531805856, previous was -13.474531805856																														
Average of 30 spectra - Confirmation Mode Complete																																
																																
<table><tr><th>Item</th><th>Min</th><th>Max</th><th>Actual</th><th>Status</th></tr><tr><td>Mass (Da)</td><td>132.8049</td><td>133.0049</td><td>132.9049</td><td>Pass</td></tr><tr><td>Sum Intensity (cps)</td><td>1000</td><td>1e+11</td><td>7558.7</td><td>Pass</td></tr></table>	Item	Min	Max	Actual	Status	Mass (Da)	132.8049	133.0049	132.9049	Pass	Sum Intensity (cps)	1000	1e+11	7558.7	Pass	<table><tr><th>Item</th><th>Min</th><th>Max</th><th>Actual</th><th>Status</th></tr><tr><td>Mass (Da)</td><td>829.4393</td><td>829.6393</td><td>829.5385</td><td>Pass</td></tr><tr><td>Sum Intensity (cps)</td><td>1000</td><td>1e+11</td><td>8221</td><td>Pass</td></tr></table>	Item	Min	Max	Actual	Status	Mass (Da)	829.4393	829.6393	829.5385	Pass	Sum Intensity (cps)	1000	1e+11	8221	Pass	
Item	Min	Max	Actual	Status																												
Mass (Da)	132.8049	133.0049	132.9049	Pass																												
Sum Intensity (cps)	1000	1e+11	7558.7	Pass																												
Item	Min	Max	Actual	Status																												
Mass (Da)	829.4393	829.6393	829.5385	Pass																												
Sum Intensity (cps)	1000	1e+11	8221	Pass																												
																																

Channel Alignment

Pass

Sum of 10 spectra

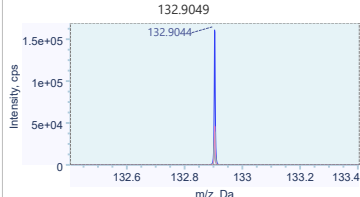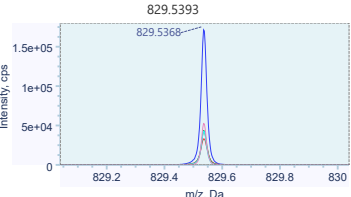									
Item	Min	Max	Actual	Status	Item	Min	Max	Actual	Status
Mass (Da)	132.6549	133.1549	132.9044	Pass	Mass (Da)	829.2893	829.7893	829.5368	Pass
Resolution	0	100000	35084	Pass	Resolution	0	100000	38241	Pass
Sum Intensity (cps)	0	1e+11	6.7117e+05	Pass	Sum Intensity (cps)	0	1e+11	1.6172e+06	Pass
Align to channel: 2 Result: Pass									
Item	Min	Max	Initial	Actual	Status				
Channel 1 vs. 2 Slope	-50	50	-0.121951398	-0.121110824	Pass				
Channel 1 vs. 2 Intercept	-50	50	-11.124393935	-10.931871514	Pass				
Channel 3 vs. 2 Slope	-50	50	0.078698755	0.071687157	Pass				
Channel 3 vs. 2 Intercept	-50	50	-4.339210041	-3.355187482	Pass				
Channel 4 vs. 2 Slope	-50	50	0.120717384	0.116474313	Pass				
Channel 4 vs. 2 Intercept	-50	50	8.491284376	9.103353985	Pass				

Instrument Name: ZenoTOF™ 7600 System
Instrument Model: ZenoTOF 7600 System
Ion Source: OptiFlow 50-200μL Micro/MicroCal
Serial Number: FB24962403


Manufacturer: AB Sciex
Target Instrument: ZenoTOF 7600 System
Optimization Set: Positive Quick Status Check

TOF MS Mass Check

Pass

Slope: 0.000490188854490918, previous was 0.000490188612503516
Delay: -13.4541511708638, previous was -13.474531805856

Average of 10 spectra - Confirmation Mode Complete

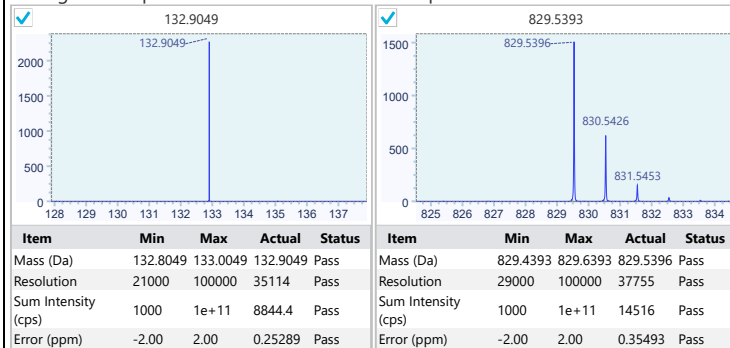


TOF MS/MS Mass Check

Pass

Slope: 0.000490188148146496, previous was 0.000490187987220103
Delay: -13.4669257821815, previous was -13.4957676613684

Average of 10 spectra - Confirmation Mode Complete

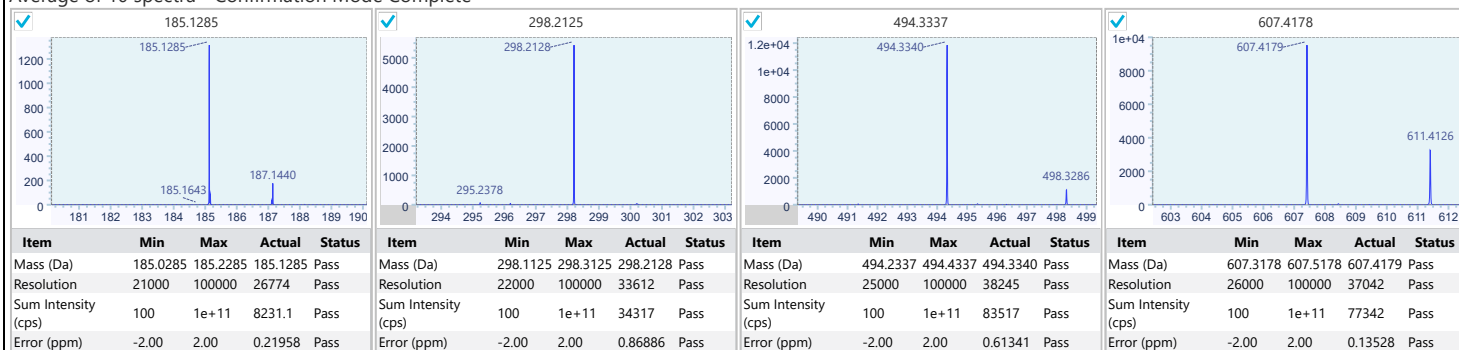

TOF MS/MS Mass Check Zeno On

Pass

Slope: 0.000490189277734504, previous was 0.000490189519273774
Delay: -13.4454189474299, previous was -13.4619974363169

Average of 10 spectra - Confirmation Mode Complete

