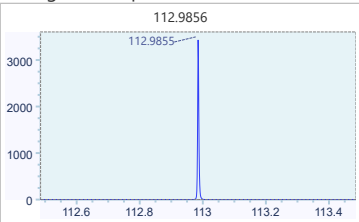
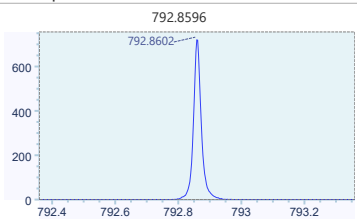
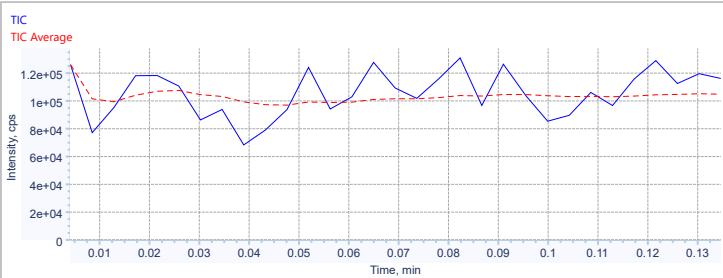


Instrument Name: ZenoTOF™ 7600 System
Instrument Model: ZenoTOF 7600 System
Ion Source: OptiFlow 50-200µL Micro/MicroCal
Serial Number: FB24962403

Manufacturer: AB Sciex
Target Instrument: ZenoTOF 7600 System
Optimization Set: Negative Quick Status Check

1. Negative MS Check

Step	Result	Details																														
Achieve Stable Spray / Modify	Pass	Slope: 0.000490038964511276, previous was 0.000490038964511276 Delay: -13.6216960433637, previous was -13.6216960433637																														
Average of 30 spectra - Confirmation Mode Complete																																
																																
<table><tr><th>Item</th><th>Min</th><th>Max</th><th>Actual</th><th>Status</th></tr><tr><td>Mass (Da)</td><td>112.8856</td><td>113.0856</td><td>112.9855</td><td>Pass</td></tr><tr><td>Sum Intensity (cps)</td><td>1000</td><td>1e+11</td><td>15432</td><td>Pass</td></tr></table>	Item	Min	Max	Actual	Status	Mass (Da)	112.8856	113.0856	112.9855	Pass	Sum Intensity (cps)	1000	1e+11	15432	Pass	<table><tr><th>Item</th><th>Min</th><th>Max</th><th>Actual</th><th>Status</th></tr><tr><td>Mass (Da)</td><td>792.7596</td><td>792.9596</td><td>792.8602</td><td>Pass</td></tr><tr><td>Sum Intensity (cps)</td><td>1000</td><td>1e+11</td><td>7597.6</td><td>Pass</td></tr></table>		Item	Min	Max	Actual	Status	Mass (Da)	792.7596	792.9596	792.8602	Pass	Sum Intensity (cps)	1000	1e+11	7597.6	Pass
Item	Min	Max	Actual	Status																												
Mass (Da)	112.8856	113.0856	112.9855	Pass																												
Sum Intensity (cps)	1000	1e+11	15432	Pass																												
Item	Min	Max	Actual	Status																												
Mass (Da)	792.7596	792.9596	792.8602	Pass																												
Sum Intensity (cps)	1000	1e+11	7597.6	Pass																												
																																

Channel Alignment

Pass

Sum of 10 spectra

Item	Min	Max	Actual	Status
Mass (Da)	112.7356	113.2356	112.9851	Pass
Resolution	0	100000	29413	Pass
Sum Intensity (cps)	0	1e+11	1.1156e+06	Pass

Item	Min	Max	Actual	Status
Mass (Da)	792.6096	793.1096	792.8583	Pass
Resolution	0	100000	37725	Pass
Sum Intensity (cps)	0	1e+11	1.1739e+06	Pass

Align to channel: 2 Result: Pass

Item	Min	Max	Initial	Actual	Status
Channel 1 vs. 2 Slope	-50	50	0.045603782	0.046240991	Pass
Channel 1 vs. 2 Intercept	-50	50	-10.044191671	-10.177094427	Pass
Channel 3 vs. 2 Slope	-50	50	-0.115200258	-0.114308094	Pass
Channel 3 vs. 2 Intercept	-50	50	-3.600434967	-3.883332087	Pass
Channel 4 vs. 2 Slope	-50	50	-0.283882493	-0.284818596	Pass
Channel 4 vs. 2 Intercept	-50	50	7.839528443	7.801046027	Pass

Instrument Name: ZenoTOF™ 7600 System
Instrument Model: ZenoTOF 7600 System
Ion Source: OptiFlow 50-200μL Micro/MicroCal
Serial Number: FB24962403


Manufacturer: AB Sciex
Target Instrument: ZenoTOF 7600 System
Optimization Set: Negative Quick Status Check

TOF MS Mass Check

Pass

Slope: 0.000490038437386946, previous was 0.000490038964511276
Delay: -13.6135203517637, previous was -13.6216960433637

Average of 10 spectra - Confirmation Mode Complete

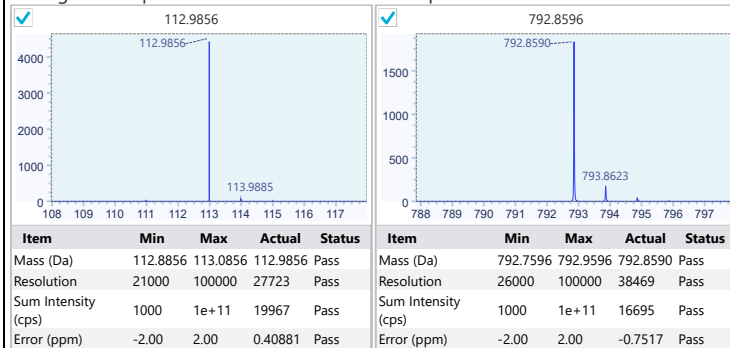


TOF MS/MS Mass Check

Pass

Slope: 0.000490037463253154, previous was 0.000490037830651855
Delay: -13.6770172812114, previous was -13.6943415111003

Average of 10 spectra - Confirmation Mode Complete

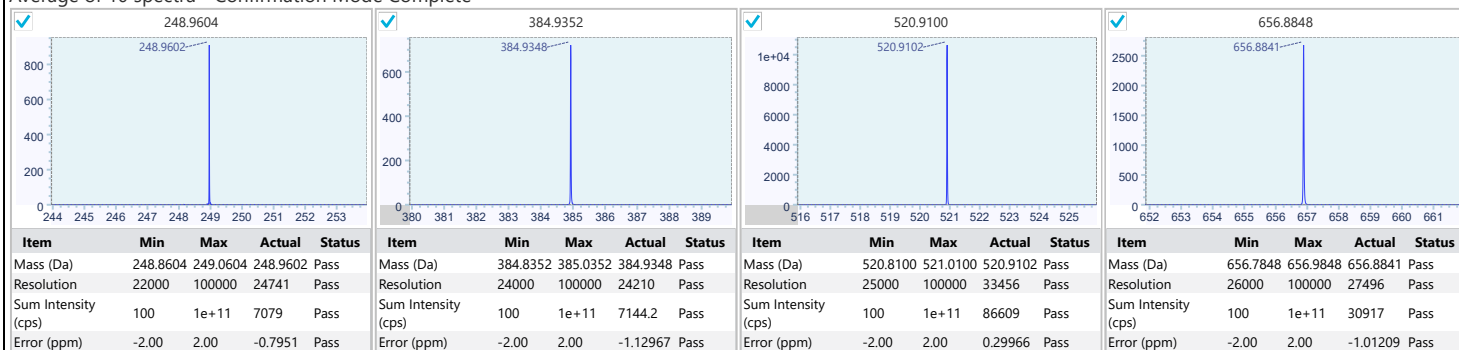

TOF MS/MS Mass Check Zeno On

Pass

Slope: 0.000490037331338447, previous was 0.00049003754350224
Delay: -13.6470264237913, previous was -13.653043613601

Average of 10 spectra - Confirmation Mode Complete

